

ADJUSTALBE
& LOCKING
FIXATION



TECHNICAL SHEET

 **comete**[®]
Control

 **comete**[®]
Extraplate



Locking system

Loop diameter: 1.35 mm*

Pull out: 1844N**



comete[®]
Control

Traction Wire :
white and green

Locking the thread splice

Loop
Thread Diameter: 1.35 mm*
Minimum length: 15 mm

Adjustment thread

Locking thread

Titanium plate TA6VELI

Length: 12 mm
Width: 4 mm
Thickness: 1.5 mm

comete[®]
Extraplate



Titanium plate TA6VELI
Length: 20 mm
Width: 7 mm
Thickness: 1.5 mm


comete[®]
Control



comete[®]
Extraplate


* Diameter measured under load

**Result obtained after dynamic test (device only) carried out as part of the CE marking of the product

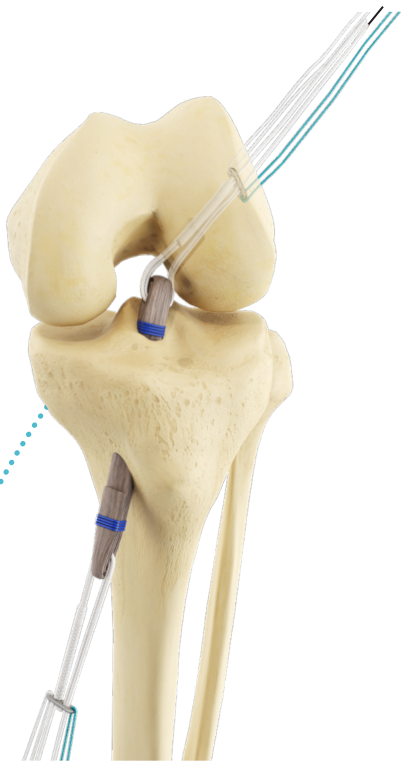

1 Femoral loop adjustment

Loop adjustment :
Alternately pulling the adjustment and locking threads
Do not force the locking thread.

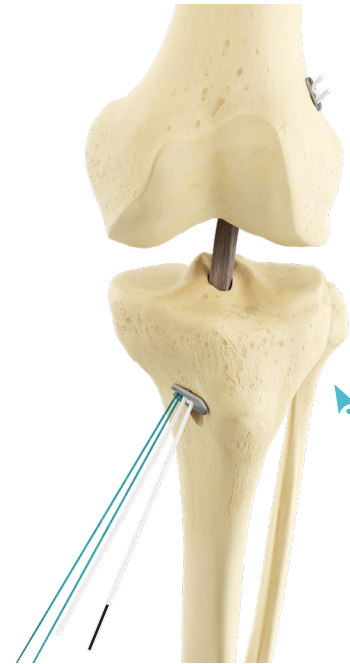

4 Tibial loop locking



2 Femoral loop locking

Locking the fixation :
Hold tension on the white adjustable thread and pull tightly the black adjustable thread to lock the loop
This action will grip the thread splice



3 Tibial loop adjustment
Graft tension adjustment

COMETE®:

comete® Control : Adjustable ACL Cortical Attachement


OAMGEFCR1U

comete® Extra plate : Plate For 6/10 mm Tunnel


OAMGEPLA1U

Amplitude instrumentation set


Modular Inside-out instrumentation set

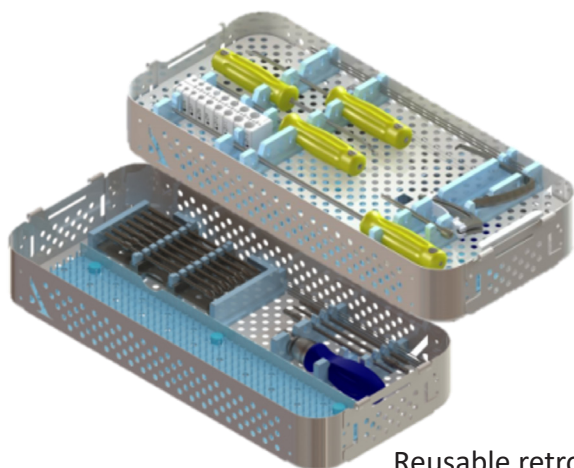

Modular Outside-in instrumentation

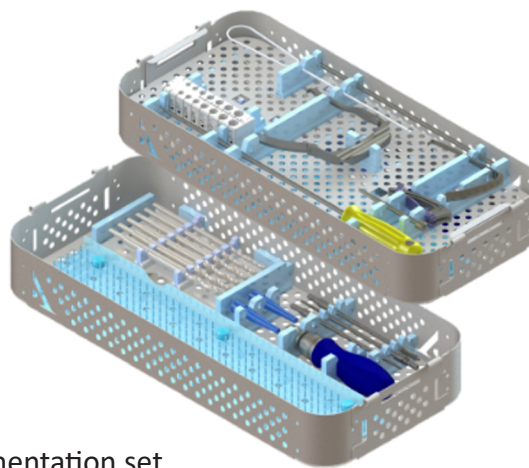

Reusable retro-drill instrumentation set

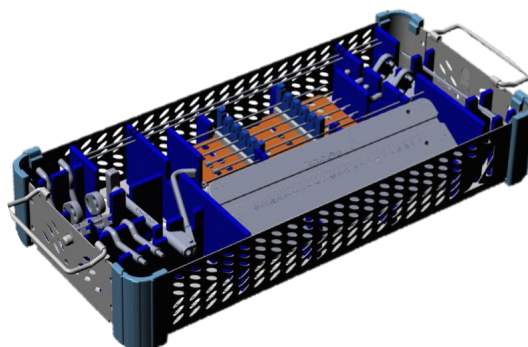

Legal manufacturer: COUSIN BIOTECH: s.a.s 8, rue de l'Abbé Bonpain, 59117 Wervicq-Sud, France

Service Clients – France :

Porte du Grand Lyon,
01700 Neyron – France
Tél. : **+33 (0)4 37 85 19 19**
Fax : +33 (0)4 37 85 19 18
E-mail : amplitude@amplitude-ortho.com

Customer Service – Export :

11, cours Jacques
Offenbach, ZA Mozart 2,
26000 Valence – France
Tél. : **+33 (0)4 75 41 87 41**
Fax : +33 (0)4 75 41 87 42